




Remote Patient Monitoring

Doccla RPM

v2.122.0

Patient/Carer App - User Manual

UDI (01)50609485900116(11)20260617(17)20280617(10)RPM21200

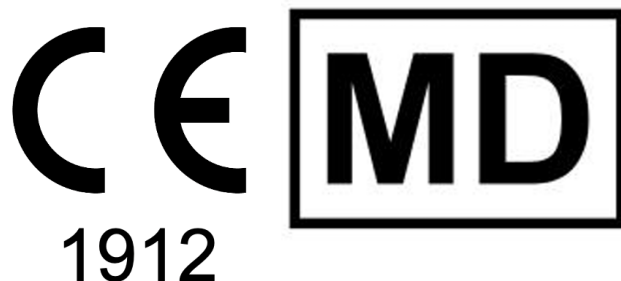

REF Doccla RPM Version: 2.122.0
Component Version: 1.43.0

Version issue date: 2026-06-17

Doccla RPM Patient App - User Manual

Author(s): Tom Austen, Rob Leyland

Table of contents

Table of contents	2
Introduction	3
Your health	4
Intended purpose	4
Medical indications	4
Measurement types in scope	4
Clinical areas in scope	5
Main functionalities	6
Benefit claims	6
Patient population	6
Intended users	7
Contraindications	7
Terms and Definitions	8
Warning and Precautions	10
Warnings	10
Precautions	10
Telemedicine platform Doccla RPM	11
Example kit	11
Supported measurements	11
Bluetooth integrated devices	11
Software integrated devices	13
App parts	14
Log on to app	14
Language selection	14
Login support - "I can't log in"	15
Forgotten login or password	15
Selecting your country	15
I have a code	16
Home menu	17
Change Password	18
Set Permissions	18
In-app Bluetooth Pairing	19
App Navigation Guidance	21
My Results	22
My Results – Temperature	22
My Replies	23
My Replies – Daily questionnaire	23
Reviewed	25

Messages	25
Start questionnaire	28
Answer a questionnaire	29
Measurement input validation	30
Information	31
Calendar	32
Video conference	32
Passive Monitoring	34
Overview	34
Setting Up Passive Monitoring	34
Checking Monitoring Status	34
Data Transmission	34
Troubleshooting	34
Important Notes	34
Carer App - Patient Selection	35
Addendum	36
Reporting patient incident	36
Precautions	36
Legal Notices	36
Manufacturer	36
Manufacturer responsibility	36
End of life (EOL)	36

Introduction

This user manual has been written for patients who must use the Doccla RPM app. The user manual has been written in a way to be thorough and understandable for patients. No prior teaching or workshop participation is necessary to understand it.

Your health



Always be mindful of your own well being and health. If you experience any new symptoms or unexpected deterioration of existing symptoms you should seek medical attention outside the use of this app. You should also contact a healthcare professional if you are unsure about your symptoms, if the app does not provide expected information, or if you are unable to use the app as intended.

Do not rely solely on this app for medical decisions.

A printed copy of this user manual can be requested free of charge by contacting Doccla UK at support@doccla.com

Intended purpose

Doccla Remote Patient Monitoring (RPM) is software intended for use in the remote monitoring of clinically stable chronic or post-operative conditions where the responsible healthcare professional has determined that intermittent remote monitoring is clinically appropriate and safe, with no risk of rapid or serious deterioration in adult and paediatric patients (≥ 1 month old) for whom intermittent remote monitoring is considered clinically appropriate and safe by the responsible healthcare professional.

This software is intended to be used by qualified healthcare professionals in non-clinical and low-acuity clinical environments, including patient homes, community care settings, rehabilitation facilities, nursing homes, and outpatient follow-up contexts.

This medical device software is contraindicated for patients who are clinically unstable, require real-time physiological monitoring, are at risk of serious or rapid deterioration, or are unsuitable for remote monitoring due to functional, cognitive, or environmental limitations.

This medical device software uses measurements obtained from compatible third-party medical devices, manual entry or questionnaires, spanning chronic disease management, post-operative care, and specialist clinical pathways in order to produce visualisations of patient data, trend displays, and alerts triggered by deviations from healthcare-professional-defined thresholds.

These outputs are intended to support clinical review, prioritisation, and follow-up actions by qualified healthcare professionals within existing clinical workflows, and to inform assessment of patient status. The outputs do not provide diagnosis, prognosis, or therapeutic recommendations and do not replace clinical judgement.

Medical indications

Doccla RPM is indicated to support the remote monitoring in the following categories of measurement types ("Measurement types in scope" section below) and clinical areas ("Clinical areas in scope" section below), either through compatible third party medical devices, manual entry or questionnaires.

Measurement types in scope

Measurement Category	Specific Parameters	Primary Clinical Area(s)
Vital signs	Blood pressure, pulse/heart rate, SpO2, respiratory rate, temperature, consciousness (AVPU)	All clinical areas (see "Clinical areas in scope" section below)
Composite early warning scores	Calculated from vital sign components, e.g. NEWS2	All clinical areas (see "Clinical areas in scope" section below)

Measurement Category	Specific Parameters	Primary Clinical Area(s)
Respiratory function	Peak flow, SpO2, FEV1, FEV6, FEV1/FEV6 ratio, FEF25-75, respiratory rate, COPD prediction score, oxygen therapy hours	Respiratory
Cardiac monitoring	ECG, pulse, LVAD parameters (pump flow, speed, power, pulse index)	Cardiology
Glycaemic monitoring	Blood glucose (mmol/L and mg/dL), blood glucose weekly average	Diabetes
Body composition & weight	Weight (kg/lb), body cell mass, fat mass, phase angle, height, BMI (derived)	Heart failure, frailty, diabetes
Activity & rehabilitation	Daily steps, daily steps weekly average, sit-to-stand test, Borg rating, VAS pain score	Frailty, respiratory, cardiology
Urinalysis	Urine blood, leukocytes, glucose, nitrite, combination dipstick, volume	Renal, specialist
Obstetrics / maternal	FeMOM ECG, FeMOM CTG, FeMOM short-term variability, FeMOM baseline, oxygen therapy, maternity BP	Obstetrics
Renal / dialysis	Peritoneal dialysis, inflow bag colour/volume, outflow bag volume, ultrafiltration volume	Renal
Haematology / biochemistry	Haemoglobin, CRP (manual entry)	Specialist
Subjective patient-reported	Multiple choice, integer/decimal scales, yes/no, free text (PROMs/symptom questionnaires)	All clinical areas (see “Clinical areas in scope” section below)

Clinical areas in scope

Clinical Area	Representative Pathway Types
Cardiology	Heart failure (optimisation, HFpEF, HFrEF, telehealth), LVAD, atrial fibrillation, arrhythmia, post-ablation, post-PCI, TIA/stroke, palpitation, community cardiology
Respiratory	COPD (including COPD prediction, early supported discharge, proactive monitoring), asthma, ARI, COVID-19, pneumothorax, pulmonary rehab, oxygen therapy, paediatric bronchiolitis, respiratory general
Frailty	General frailty, care home, frailty with/without observations, frailty falls, proactive frailty
Diabetes	Type 1/2 blood glucose monitoring, gestational diabetes (diet, HBGM, insulin, metformin), diabetes oncology
Obstetrics	Hypertension screening, blood pressure management, maternity BP, preeclampsia, ROCH blood pressure
General	Step-up/step-down, home monitoring, virtual wards, community matrons, MOA early supported discharge
Palliative	Palliative care, comfort tracker (long/short), IPOS, care home
Renal	UTI, urology, peritoneal dialysis
Specialist	OPAT, paediatrics (general, respiratory, oncology, haematology), orthopaedics, surgery, colorectal, hepatology, infectious diseases, MPox, gynaecology

Main functionalities

Data Collection and Integration: Acquisition of physiological and non-physiological data through compatible third party medical devices, manual entry, or connected systems.

Alert Configuration and Generation: Allows healthcare professionals to configure thresholds and/or validated algorithms that generate alerts when patient data deviate from predefined criteria. Where configured, this includes calculation of composite early warning scores from patient-submitted vital signs data using clinician-defined scale settings. Alerts are based solely on clinician-defined inputs and do not interpret, diagnose, or recommend treatment.

Data Visualisation and Trend Analysis: Presentation of historical and current patient data to support clinical interpretation and prioritisation.

Doccla RPM supports intermittent monitoring only and is not intended for real-time monitoring.

Benefit claims

Doccla RPM is intended to deliver clinical benefit by enabling timely clinician response to deviations from clinician-defined thresholds in clinically stable patients, and by providing a monitoring experience that patients find acceptable and are willing to sustain. The two clinical benefits being formally claimed, and their corresponding outcome parameters, are set out below.

Benefit	Benefit Claim	Mechanism	Outcome Parameters
Benefit 1	Timely clinical review and prioritisation	Doccla RPM generates alerts when measurements from clinically stable patients breach clinician-defined thresholds, enabling timely review and documented follow-up action within existing clinical workflows	Time from threshold breach to first qualifying clinician review. Proportion of alerts with a documented clinical action within a prespecified time window. Time to first documented action.
Benefit 2	Acceptable and usable remote monitoring	Patients are willing and able to engage with remote monitoring over the duration of their care pathway, supporting the completeness and clinical utility of monitoring data	Patient-reported acceptability and usability. Patient engagement and retention. Compliance with expected measurement schedules.

Patient population

Adult and paediatric patients (≥ 1 month old) for whom intermittent remote monitoring is considered clinically appropriate and safe by the responsible healthcare professional.

Intended users

Doccla RPM is intended to be used by:

- Qualified healthcare professionals: For configuring monitoring parameters, reviewing patient data, and responding to alerts through the clinician portal.
- Patients and caregivers: For entering or transmitting data via the patient or carer mobile application.

All users must receive appropriate training in the safe operation of the software and any connected medical devices.

Contraindications

Doccla RPM is contraindicated for patients who are clinically unstable, require real-time physiological monitoring, are at risk of rapid or serious deterioration, or are unsuitable for remote monitoring due to functional, cognitive, or environmental limitations.

The above contraindications are further detailed below by clinical area, based on results of risk management. Suitability for remote monitoring must be assessed prior to enrolment and reviewed periodically by the responsible healthcare professional.

The contraindications below describe clinical presentations that are outside the intended purpose of Doccla RPM. Where a clinical area is named, the contraindication applies to the acute or unstable presentation of that condition. The diagnosis alone does not determine suitability - the responsible healthcare professional's assessment of the patient's current clinical status does.

Clinical Area	Contraindication
Cardiology	Patients in acute cardiac decompensation (e.g., acute decompensated heart failure, haemodynamically unstable arrhythmia, acute MI) requiring inpatient or emergency management.
	Patients requiring real-time cardiac monitoring or surveillance.
	Post-operative cardiac patients prior to clinical stabilisation where condition has not yet been confirmed stable by the responsible clinician.
Respiratory	Acute exacerbations of COPD or asthma requiring inpatient treatment or active titration of unstable oxygen requirements.
	Acute respiratory failure or patients requiring high-flow oxygen therapy where condition cannot be managed remotely.
Frailty	Patients with severe physical frailty who are unable to operate the patient application or connected measurement devices independently, where carer support is not available or cannot be confirmed.
	Patients with significant cognitive impairment or severe mental health conditions who are unable to engage with the monitoring process, where carer-assisted use has been assessed as insufficient.
	Patients assessed as high risk for rapid deterioration associated with falls or acute frailty episodes, where the intermittent monitoring cadence of Doccla RPM is insufficient to detect deterioration in time.
Diabetes	Patients with poorly controlled Type 1 diabetes requiring continuous or near-continuous glucose monitoring integrated with insulin delivery (e.g., closed-loop systems).
	Diabetic ketoacidosis (DKA) or hyperosmolar hyperglycaemic state (HHS) - acute conditions requiring immediate clinical intervention.
	Patients with gestational diabetes where the rapidity of physiological change in glycaemic status and insulin sensitivity requires direct, continuous clinical assessment at the point of care.
Obstetrics	Hypertensive crisis requiring emergency assessment.

Clinical Area	Contraindication
	Pre-eclampsia with severe features (e.g., severe-range blood pressures, symptoms of end-organ involvement) requiring inpatient management.
	Patients where blood pressure variability or instability is such that intermittent readings cannot provide a clinically safe picture of the patient's status.
General	Patients requiring monitoring at a frequency or acuity level that exceeds Doccla RPM's intermittent monitoring capability.
	Patients without access to a compatible smartphone or tablet, or without sufficient and reliable internet connectivity, where this cannot be resolved at onboarding and neither manual data entry nor carer-assisted data submission has been assessed as a clinically acceptable alternative."
	Patients who are unable or unwilling to use Doccla-supported measurement devices, and for whom manual data entry is not a clinically acceptable alternative.
	Patients who cannot meaningfully consent to remote monitoring, including understanding the non-emergency nature of the service and associated data collection.
Palliative	Patients in the terminal phase of illness where remote monitoring has not been specifically agreed with the patient, their family, and the clinical team as consistent with their goals of care.
	Patients for whom alert-based threshold monitoring is not aligned with their agreed care plan (e.g., where a do-not-escalate decision is in place but thresholds have not been adjusted accordingly).
Renal	Patients with acute kidney injury (AKI) or rapidly deteriorating renal function requiring frequent in-person clinical review.
	Patients undergoing acute or urgent-start peritoneal dialysis where fluid balance and electrolyte management requires more frequent clinical input than intermittent remote monitoring can support.
Specialist	Paediatric patients under 1 month of age (minimum age for intended use).
	Post-operative surgical or orthopaedic patients in the immediate post-operative period prior to clinical stabilisation, where wound status or acute complications require in-person assessment.
	Patients on OPAT (outpatient parenteral antibiotic therapy) or infectious disease pathways where systemic infection may deteriorate rapidly and requires a monitoring cadence beyond Doccla RPM's intermittent model.
	Oncology and haematology patients undergoing active treatment where the treatment protocol or current clinical status creates a risk of acute life-threatening complications requiring inpatient management

Terms and Definitions

Term / Abbreviation	Definition
AKI	Acute kidney injury
ARI	Acute respiratory infection
AVPU	Alert, Voice, Pain, Unresponsive - a scale used to assess level of consciousness
BMI	Body mass index
BP	Blood pressure
COPD	Chronic obstructive pulmonary disease

Term / Abbreviation	Definition
CRP	C-reactive protein - a biochemical marker of inflammation
CTG	Cardiotocography - monitoring of fetal heart rate and uterine contractions
DKA	Diabetic ketoacidosis
ECG	Electrocardiogram
FEF25-75	Forced expiratory flow between 25% and 75% of forced vital capacity
FeMOM	Fetal and Maternal Outpatient Monitoring - a proprietary obstetric monitoring device
FEV1	Forced expiratory volume in one second
FEV6	Forced expiratory volume in six seconds
FEV1/FEV6	Ratio of FEV1 to FEV6, used as a surrogate for FEV1/FVC in spirometry
FVC	Forced vital capacity
HBGM	Home blood glucose monitoring
HFpEF	Heart failure with preserved ejection fraction
HFrEF	Heart failure with reduced ejection fraction
HHS	Hyperosmolar hyperglycaemic state
HCP	Healthcare professional
IPOS	Integrated Palliative care Outcome Scale
LVAD	Left ventricular assist device
MDR	EU Medical Device Regulation 2017/745
MOA	Medicine for Older Adults
NEWS2	National Early Warning Score 2 - a standardised clinical scoring system for detecting patient deterioration based on vital signs
OPAT	Outpatient parenteral antibiotic therapy
PCI	Percutaneous coronary intervention
PROMs	Patient-reported outcome measures
ROCH	Remote optimisation of chronic hypertension
RPM	Remote patient monitoring
SaMD	Software as a Medical Device
SpO2	Peripheral oxygen saturation, measured by pulse oximetry
TIA	Transient ischaemic attack
UTI	Urinary tract infection
VAS	Visual analogue scale
WCAG	Web Content Accessibility Guidelines

Warning and Precautions

Warnings

- Doccla RPM must not be used as a substitute for emergency medical services. If a patient's condition deteriorates acutely or requires immediate clinical intervention, contact emergency services without delay.
- Doccla RPM supports intermittent monitoring only. It must not be used for real-time monitoring of patient status. Alerts are threshold-based and do not guarantee detection of all clinically significant changes.
- Do not attempt to modify, reverse-engineer, or alter any component of the Doccla RPM software. Unauthorised modification may compromise the safety and performance of the device and renders any CE marking invalid.
- Doccla RPM must only be accessed using login credentials assigned to the individual user. Credentials must not be shared with or disclosed to any other person

Precautions

- Ensure that Bluetooth is enabled on the patient's device during the monitoring period. Where Bluetooth connectivity cannot be established, manual data entry or carer-assisted data submission should be considered as alternatives. Patients and carers must be instructed on all available data submission methods at onboarding.
- All users must receive appropriate training before using Doccla RPM. Clinicians must be trained by a qualified representative or a colleague who has completed such training. Patients and carers must be onboarded by a trained member of the clinical or support team.
- Only measurement devices listed as compatible with Doccla RPM must be used for automated data capture. Use of unsupported devices may result in data not being received or being incorrectly parsed. Where automated capture is not possible, manual data entry must be assessed for clinical acceptability.
- Patients and carers must be clearly informed at onboarding that Doccla RPM is a monitoring support tool and is not an emergency service. Patients must understand that the absence of an alert does not confirm that their condition is stable or that no clinical action is required.
- Endpoint devices used to access Doccla RPM - including clinician workstations, patient smartphones, and carer devices - must be maintained with current operating system updates and up-to-date endpoint security software.

Telemedicine platform Doccla RPM

In Doccla RPM you will be using a mobile device (tablet or phone) along with an optional set of medical measurement devices: The combination of mobile device + measurement device will be referred to as an Doccla RPM kit in this user manual.

Example kit

The example kit consists of the following (yours may vary):

- A tablet, with the Doccla RPM Patient App installed,
- a blood pressure monitor,
- a thermometer,
- a pulse oximeter.



Supported measurements

The Doccla RPM app supports manual input of measurements as well as automatic input of measurements using Bluetooth enabled devices. The following set of devices are all supported by Doccla RPM:

Bluetooth integrated devices

Device type	Manufacturer	Model name	Device Image
Activity Tracker	Beurer	AS80	

Blood pressure monitor	A&D	UA-651BLE	
Blood pressure monitor	A&D	UA-656BLE	
Blood pressure monitor	TeleRPM	BB2284-AE01	
ECG	AliveCor	KardiaMobile6L	
Multi-Vital ECG Sensor	Vivalink	VV350	
Multi-Vital Bracelet	Corsano	CardioWatch 287-2B	
Oximeter	Masimo	MightySat Rx	
Oximeter	Nonin	3150BLE	
Oximeter	Nonin	Onyx 3230	
Oximeter	Kinetik	MD300CI218R	

Peak Flow Meter	Vitalograph	asma-1 BT	
Thermometer	Kinetik	ET-100EB	
Weight scale	A&D	UC-352BLE	
Weight scale	Transtek	Body Scale GBS-2012-B	

Software integrated devices

Measurement type	Manufacturer	Model name
Smartphone Camera PPG	Happitech	Heart Rhythm SDK


App parts

The following section describes different parts of the Doccla Patient App.

Log on to app

When you open the Doccla app you will be presented with the login screen, enter your username and password to log in. If you check the “Remember me” checkbox, the app will remember your password and you will not need to enter your username and password again unless you explicitly log out from the app.

If using the “Remember me” function please keep in mind that personally sensitive data is visible through the app and do not share your username and password with anyone.

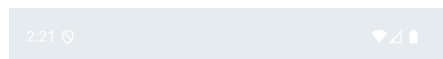

doccla


Log in to continue



Remember me ☒


Log in

By logging in you agree to our [privacy policy](#)

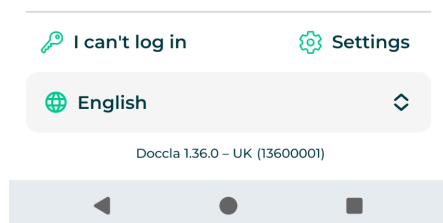

Language selection


There is a language selector at the bottom of the login screen. To change the language press onto the language indicator that opens a menu of supported languages, press onto a language to select it.

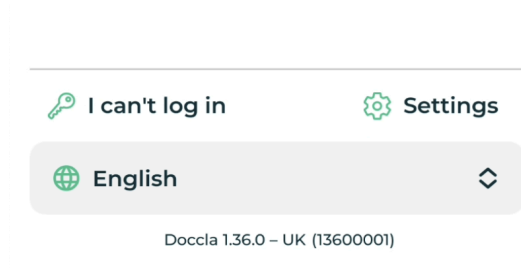

Login support - “I can’t log in”


There is an option for users that are struggling to log in to the app below the “Privacy Policy”.

When pressed the user is presented with options.

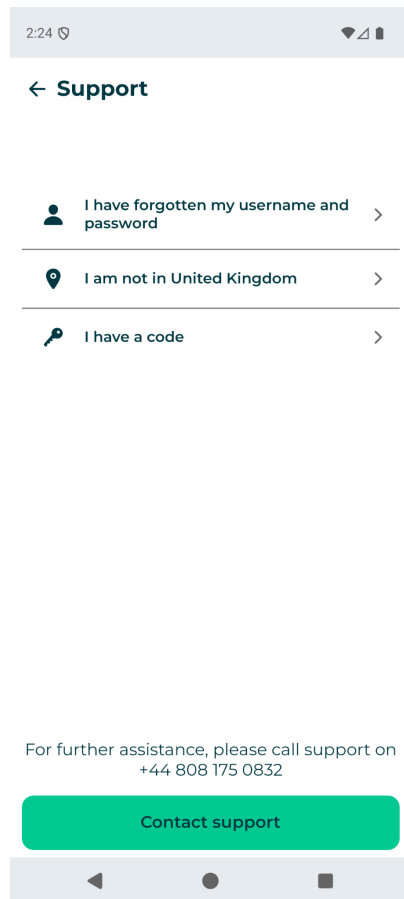


Forgotten login or password

When “I have forgotten my username and password” is pressed the instructions to call Customer Support are displayed.

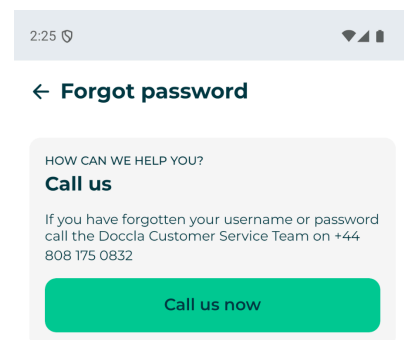



Selecting your country

When “I am not in [country name]” is pressed the user can select which country they are based in and press “Save Location”. When “Save Location” is pressed the user is returned to the Login Screen.

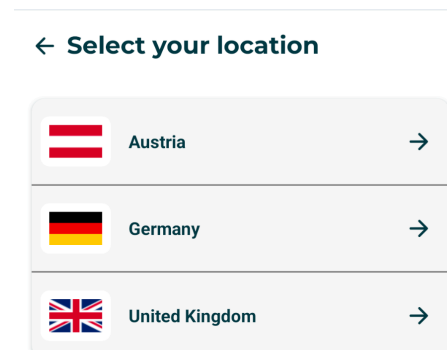

I have a code

If you have been given a QR code you can scan this and it will direct you to the login page.

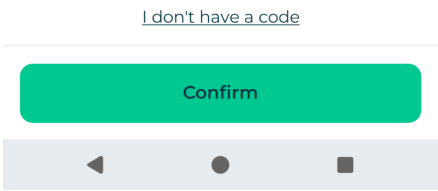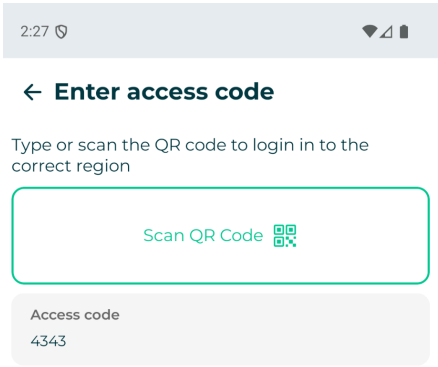

Home menu

From the home menu you can navigate to one of the following submenus by pressing their corresponding button:

Start questionnaire: Press to show all questionnaires that can or need to be completed.

Start ECG Measurement: Press to open the ECG submission module.

Messages: Press to show messages between you and the healthcare professionals.

Reviewed: Press to show acknowledgements of your recent measurements by healthcare professionals.

My results: Press to show tables and graphs of measurement results previously submitted by you.


My replies: Press to show tables of questionnaire replies previously submitted by you.

Information: Press to show links to further relevant material for your condition on the web.


Calendar: See upcoming video conferences between you and a healthcare professional.

The menu items that have a blue background and a number next to the title, here **Start questionnaire (2)** and **Messages (3)**, indicate that there are questionnaires to be answered and new messages from the clinicians, respectively.

Note: Depending on your condition you may not necessarily see all these menu items when logging into the Doccla app.

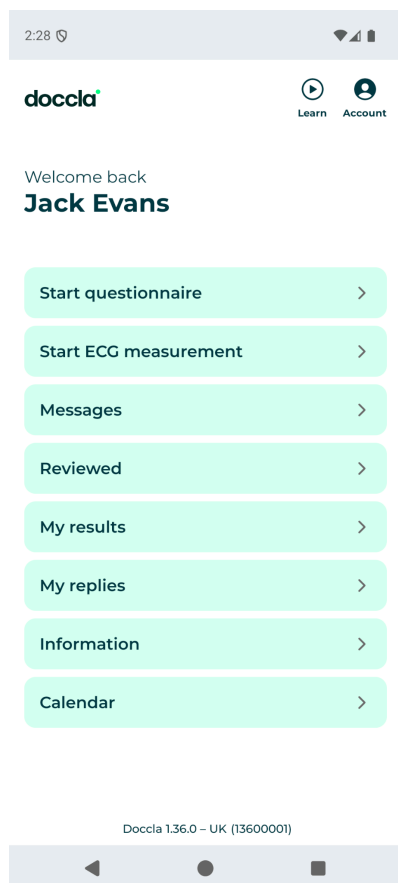

Press the **Account** button  in the top right hand corner to access the app **Settings**.

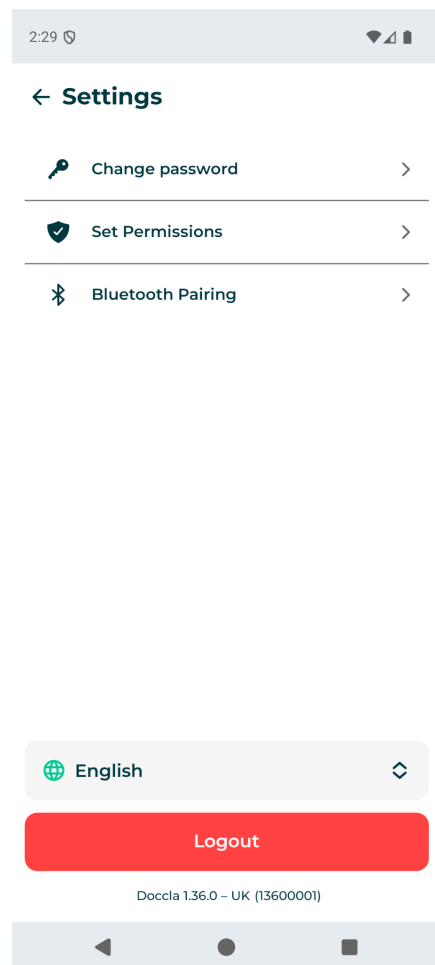

- Press the **Change password** button to go to the change password page.
- Press the **Set Permissions** button to manage the permissions granted to the app.
- Press the **Bluetooth Pairing** button to pair a medical device from the supported device list.
- Press the **Select Language** button to change the language of the app.
- Press the **Logout** button to log out of the app.

Change Password


This menu allows you to change your password by entering your current password followed by a new password twice and pressing change.

You can always press the **back** button to return to the previous menu.

Set Permissions

All of the permissions the app requires are listed shown with a green tick  if the permission has been granted or with a red cross  if the permission has not been granted.

To enable a permission click on the red cross  and on a pop up you can select what level of permission to grant.

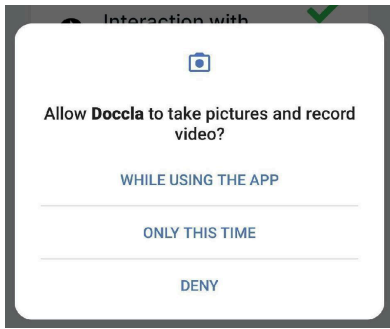


In-app Bluetooth Pairing

The Bluetooth pairing option **only** allows you to pair medical devices integrated with the Doccla app, these devices are listed on pages 11 to 13 of this User Manual.

After selecting the Bluetooth pairing option the app will display any devices that are already paired with the smartphone or tablet and offer the option to “Begin Scanning” for devices that are in pairing mode.


← Bluetooth Pairing

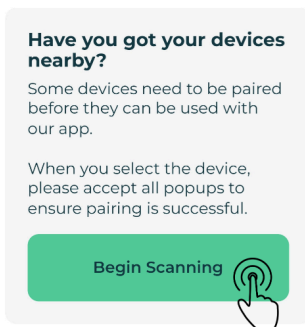

Paired devices

No paired devices found

When pressing “Begin Scanning” all nearby devices that are in pairing mode (and not already paired with the smartphone or tablet) will appear. To pair with that device select it from the list by pressing the **+** button next to the device’s name and follow the on-screen instructions.

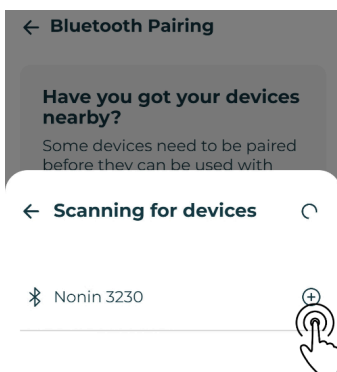

On Android smartphones and tablets you will see that the device now appears in the “Paired devices” list. On Apple smartphones and tablets (running iOS) you will need to navigate to the smartphone/tablet’s Settings app to review paired Bluetooth devices.

Have you got your devices nearby?

Some devices need to be paired before they can be used with our app.

When you select the device, please accept all popups to ensure pairing is successful.

Begin Scanning

Paired devices


 Nonin 3230 

App Navigation Guidance



Video guidance is available for how to use application features by pressing on the Learn icon **Learn** at the top of the screen, for example when inputting a blood pressure measurement in a questionnaire.

When pressed the icon will open a video window with a red play  button, when this is pressed the video will start.

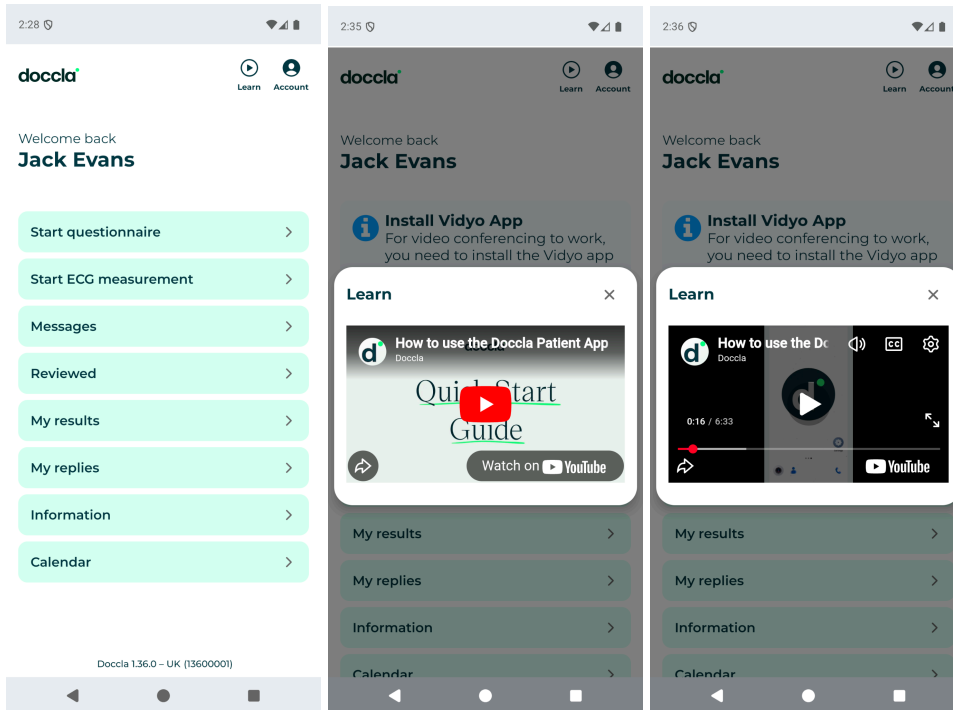


When the video has started it can be

- Made full screen by pressing the square button in the bottom right corner of the video 
- Paused by pressing the two lines button in the bottom left corner of the video 
- Closed by pressing the X button to the top right of the video window 

My Results

This menu allows you to see previous submitted results and measurements as tables and graphs. Press one of the measurement types in the list to see more details.

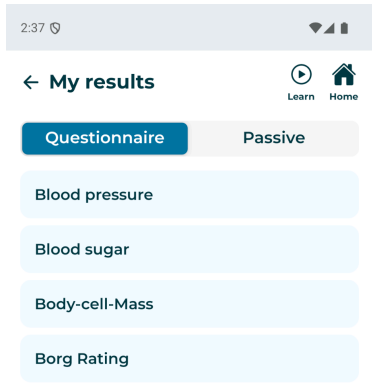


You can always press the **home** button to return to the home menu.

My Results – Temperature

This is an example of results for a specific measurement type that was chosen in the previous **My results** menu. Here you can see a graph and a table of results, where you can filter the amount of results shown by pressing on one of the time filters at the top of the menu.


You can always press the **back** button to return to the previous menu.

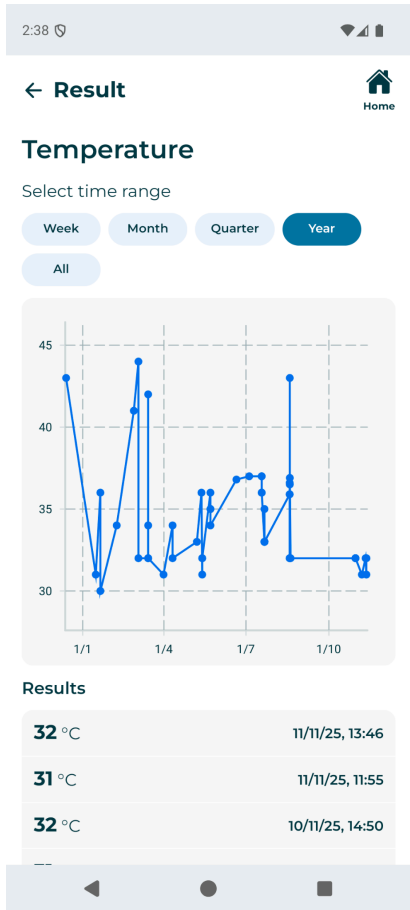

My Replies

This menu allows you to see previous submitted questionnaire replies displayed as tables. Press one of the questionnaire replies in the list to see more details.

You can always press the **home** button to return to the home menu.

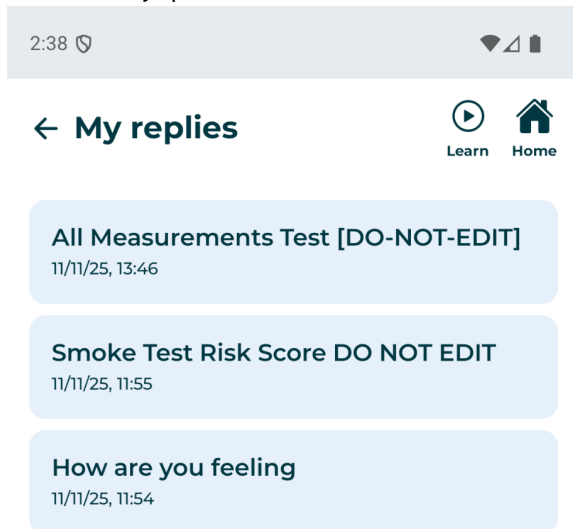

My Replies – Daily questionnaire

This is an example of a reply for a questionnaire that was chosen in the previous **My replies** menu. Here you can see a table with a column for the **Questions** asked and the **Replies** given.

You can always press the **back** button to return to the previous menu.

[← Reply](#)

YOUR ANSWERS

Daily

Submitted at 20/08/25, 12:06

Please enter your step count for today

Daily steps: 56 step count

Please enter your pulse below

Pulse: 56 BPM

Have you experienced unusual or uncomfortable heartbeats or palpitations in the past 24 hours?

Yes

How many episodes of unusual or uncomfortable heartbeats or palpitations have you experienced in the past 24 hours?

1 to 5 episodes

Did you feel anxious with these episodes of unusual or uncomfortable heartbeats or palpitations?

Occasionally

Unusual Heartbeats - AHA Questionnaire

3

Have you been short of breath during normal activities in the past 24 hours?


Yes

Reviewed

This menu shows you when a healthcare professional has reviewed your most recent results. A list of messages is displayed, one for each of your recent results sent in to the healthcare professionals including the date and time at which it has been reviewed.

If there are no results reviewed yet, you will get the message “There are currently no reviewed results to show”.

You can always press the **home** button to return to the home menu.

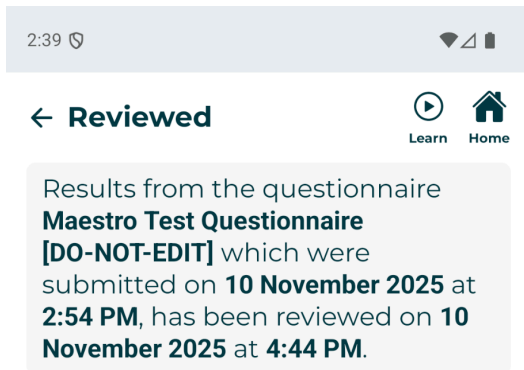


Messages

This menu shows the messages sent between you and the healthcare professionals. The newest message will be at the bottom of the screen.

At the bottom of the screen in the white input field, you can enter a message for the healthcare professionals.


After writing the message press the  button to send it to the healthcare professionals. You can always

press the **back** button  to return to the previous menu.


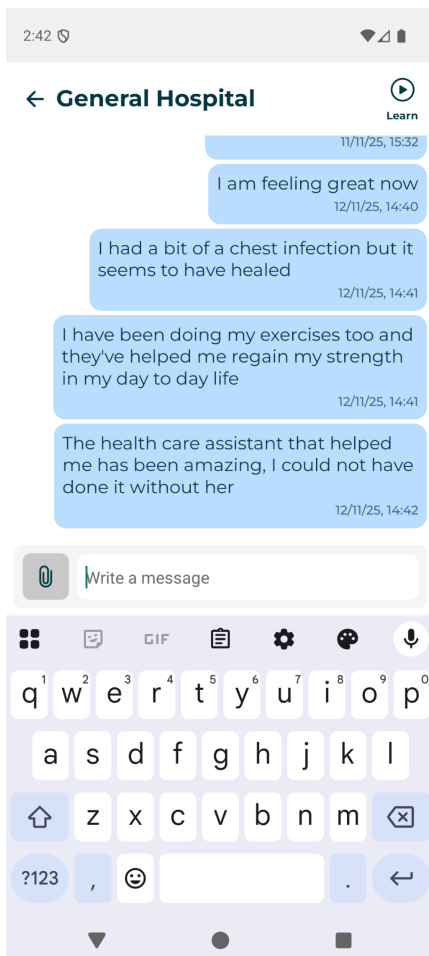




To attach an image to a message, click the button next to the send button. Depending on the mobile device, you can choose to take an image with the device's camera, or upload an existing one.

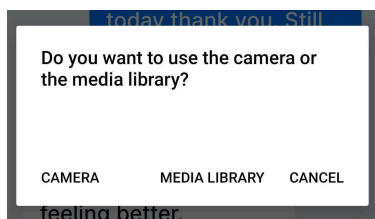

After choosing an image, a preview will be shown next to the message.

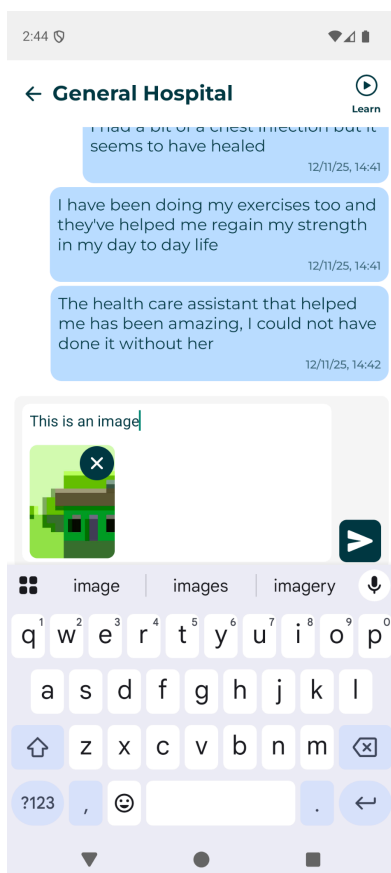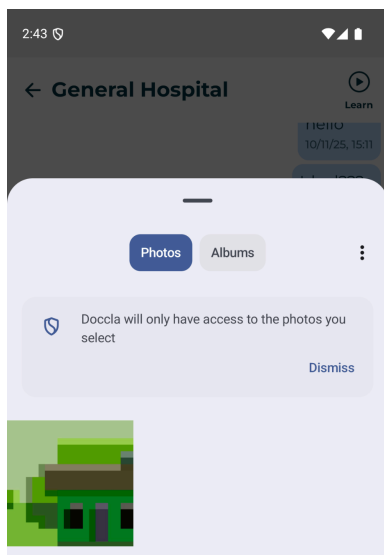

When a sent message contains an image, a thumbnail will be shown in the conversation. Click to enlarge.




If an image is deleted, a gray placeholder will be shown instead of the original image.


Start questionnaire

This menu shows all the questionnaires currently available for you to answer. Press one of the questionnaires in order to start completing it.


You can always press the **home** button to return to the home menu.

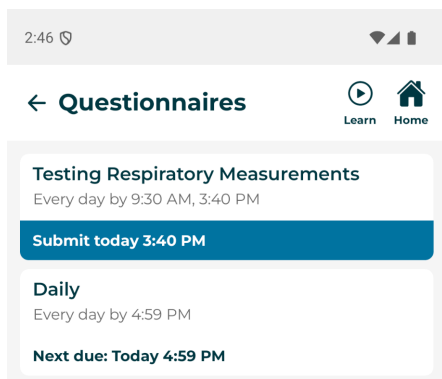


The frequency you need to submit your questionnaire will be displayed under the name of the questionnaire.

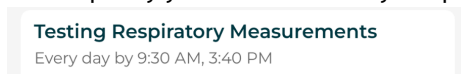

The questionnaires with a blue row indicate that these should be answered before the deadline displayed underneath the questionnaire title, e.g. the **Testing Respiratory Measurements** questionnaire has a deadline at *Today, 3:40PM*.

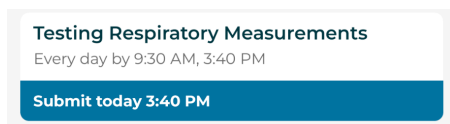

If a questionnaire has a visible deadline but has not been marked with a blue background yet, that means you should not complete the questionnaire just yet but wait until you get closer to the deadline.

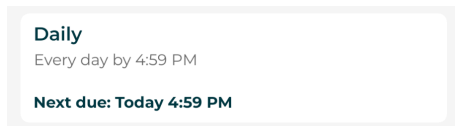

Answer a questionnaire

This is an example of a questionnaire that was chosen in the previous screen.

When answering a questionnaire, you are led through a series of different questions, each either showing a text, a simple input field, a measurement input field or a multiple choice question.

For each question, you navigate to the next by clicking on one of the buttons at the bottom of the screen.

The number of buttons may vary between the different screens, i.e. a **Next** button, a **No** and **Yes** button or a **Skip** and **Next** button.

At the end of a questionnaire you must choose to either send the questionnaire results to the healthcare professionals or omit the results.

When asked to perform a measurement please find the appropriate measurement device and take a measurement. The medical devices should already have been paired with the tablet via Bluetooth allowing for automatic measurements.

If you are using devices that are not using Bluetooth, you enter the result into the questionnaire answer box.

You can always press the **back** button to go back to the previous question to choose a different answer or the **home** button to return to the home menu.



Enter the value in the input field and press **Next**. (Consult the help text if needed, shown when clicking on the **Help** at the top).

To send the completed questionnaire results to the healthcare professionals, press **Send answers**.


An acknowledgment is shown when the completed questionnaire is successfully received.

The image displays three sequential screens of a mobile application used for answering a questionnaire.

- Screen 1:** Features a back arrow and a home icon at the top. The main text reads "IMPORTANT" followed by "Please gather the medical equipment provided and make sure you are sat calmly in an upright position". A blue "Next" button is at the bottom.
- Screen 2:** Features a back arrow and a home icon at the top. The main text reads "Please take your oxygen saturation & pulse with the oximeter provided and record your measurements below". Below this are two input fields: "Saturation" with the value "98" and "Pulse" with the value "62". A blue "Next" button is at the bottom.
- Screen 3:** Features a back arrow and a home icon at the top. The main text reads "All questions answered.". A green "Send answers" button is at the bottom.


Measurement input validation

When answering a questionnaire there is validation on the values that can be entered for measurements. If a value is input outside of the expected limits for a measurement a warning message will appear with the upper and lower limits expected for that measurement type.

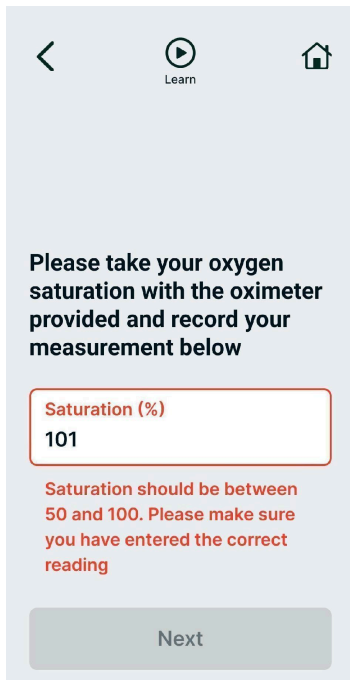

The screenshot shows a mobile app interface for recording oxygen saturation. At the top, there are navigation icons: a back arrow, a play button labeled 'Learn', and a home icon. Below these is a light blue header area with the text: "Please take your oxygen saturation with the oximeter provided and record your measurement below". A text input field is labeled "Saturation (%)" and contains the value "101". Below the input field, a red warning message reads: "Saturation should be between 50 and 100. Please make sure you have entered the correct reading". At the bottom, there is a greyed-out button labeled "Next".


The limits are designed to screen unlikely values, so the “Next” button will be greyed out and disabled until the value is either changed or three seconds have elapsed when the “Next” button turns blue and can be pressed to progress.

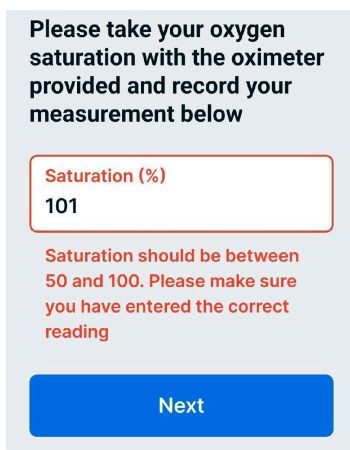

This screenshot shows the same app interface as the previous one, but the "Next" button at the bottom is now blue and active. The warning message and the input field value remain the same.

Information


This menu shows a list of links to pages, videos or documents containing further information that might help you in your treatment.

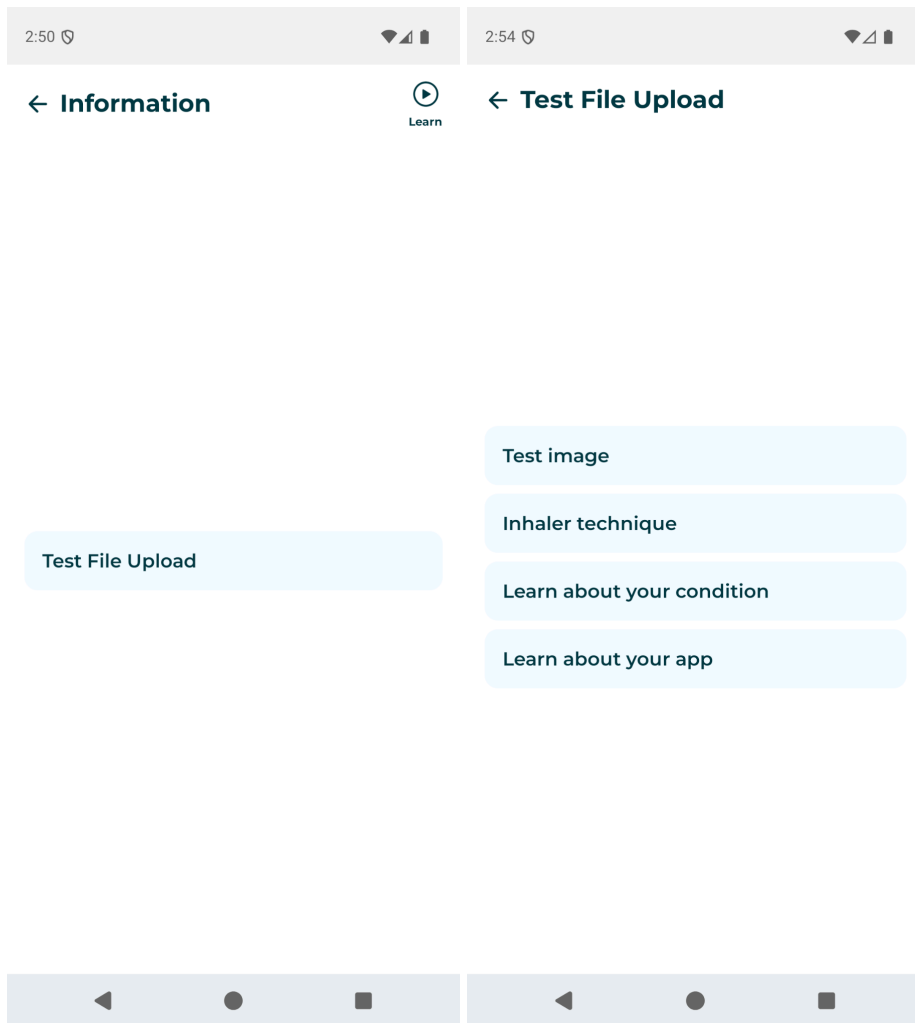

When clicking such a link be aware that you are leaving the Doccla app and opening the link in the web browser on your mobile device.


Calendar

The menu allows you to view upcoming and past video conferences with a healthcare professional.

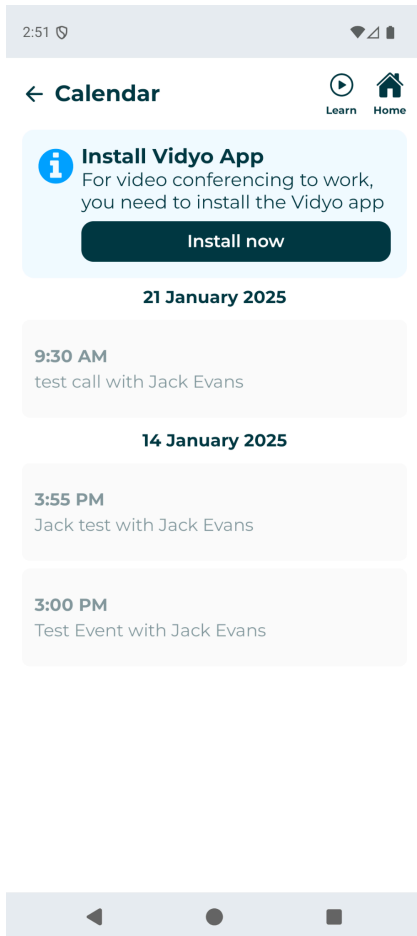

Video conference


Your clinician can choose to start a video conference with you. They will call you up on your mobile device, after which you have to accept the call by pressing the **Join video conference** button.

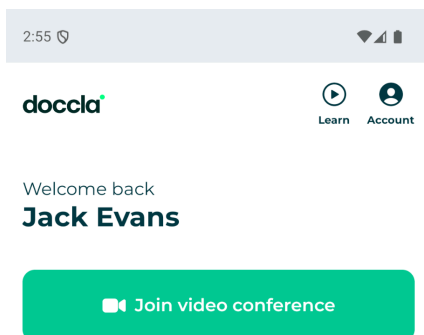

- The video conference will open in the VidyoConnect app and you will need to click **Join**.
- You will see the clinician who has initiated the call in the centre of the screen and in the top right hand corner you will see a preview of yourself. You use the controls below
 - You can turn your camera on/off by pressing the button 
 - You can turn your microphone on/off by pressing the button 
 - You can end the call by pressing the button 
 - **Do not use** the functionality of the button  is not supported.

- When the call has ended you can return to the Doccla Patient App by clicking  at the bottom of the Doccla RPM Patient App - User Manual v1.43 - v2.122.0

screen.

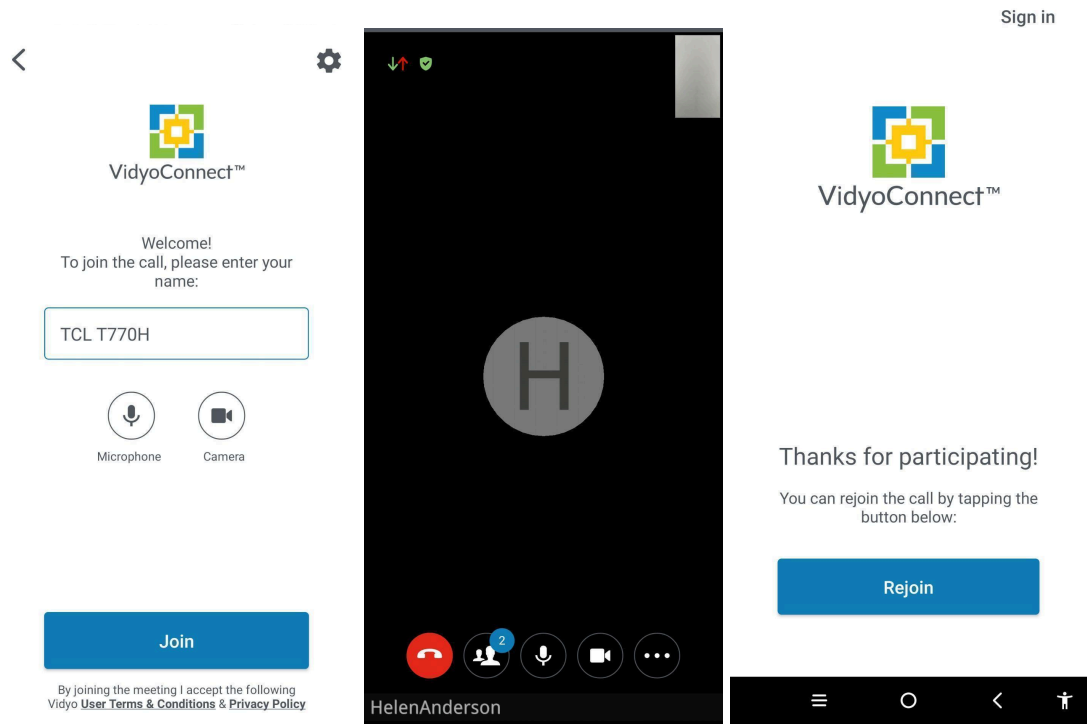

Passive Monitoring

Overview

Passive Monitoring allows the Doccla app to automatically collect physiological data from a compatible wearable device without requiring manual input. This feature enhances patient monitoring by passively tracking key vital signs and transmitting data to healthcare professionals.

Setting Up Passive Monitoring

1. Prepare the wearable device:
 - Ensure the wearable is charged. Plug the provided charger into a power adapter and charge the device.
2. Pairing the wearable:
 - Open the Doccla App and navigate to the home screen.
 - Press “Pair”
 - The app will search for the wearable. A message will appear:
 - i. Now remove the wearable from the charger
 - ii. “Bracelet found, start pairing process.”
 - If the wearable is not found, ensure Bluetooth is enabled on your mobile device.
 - iii. Once found, follow on-screen instructions to complete pairing.
3. Complete pairing profile.
4. Complete wearing optimisation
 - Follow on-screen instructions
 - i. Remain still and rest your arm on a flat surface while completing this step.
5. Monitoring Activation:
 - Once paired, Passive Monitoring starts automatically.
 - The app will notify you if data is not being recorded (e.g., if the bracelet is disconnected or battery is low).

Checking Monitoring Status

- The wearable device status is displayed on the Home Screen.
 - If monitoring is active, you will see a status indicating that data is being recorded.
 - If monitoring is not active, the app will display a “Not recording!” message.
- To navigate to the wearable page in the app press on the wearable card.
- Ensure the bracelet is worn correctly.

Data Transmission

- Data is automatically sent to healthcare professionals.
- No manual data entry is required.
- Healthcare professionals can review trends over time.

Troubleshooting


- If pairing fails:
 - Make sure the bracelet is charged.
 - Restart the bracelet by placing it back on the charger for a few seconds.
- If monitoring stops:
 - Ensure the bracelet is snug on your wrist.
 - Reconnect in the wearable device section.

Important Notes

- Passive Monitoring is not a substitute for emergency care. If you experience symptoms, seek medical assistance.
- The bracelet does not provide real-time alerts for critical conditions.
- Always ensure the bracelet is charged and correctly positioned on your wrist.
- If the device is not on the latest firmware (“Update available” will be displayed) then press to update.

Carer App - Patient Selection

When you have logged in you will be presented with a list of the patients that you care for.

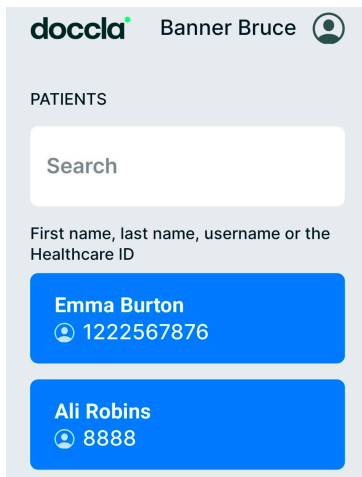

doccla Banner Bruce

PATIENTS

Search


First name, last name, username or the Healthcare ID

Emma Burton
1222567876

Ali Robins
8888

It is recommended that you search for the patient in the search box at the top of the Home Screen. You can search for patients with any of the below fields from the Clinician Dashboard -

- First name
- Last name
- Healthcare ID
- Username

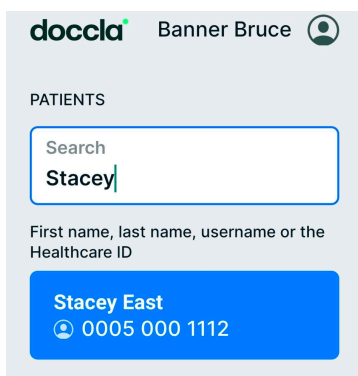

doccla Banner Bruce


PATIENTS

Search
Stacey

First name, last name, username or the Healthcare ID

Stacey East
0005 000 1112

You will be asked to confirm that you have selected the correct patient, it is important to check the Healthcare ID in the app matches the Healthcare ID on the Clinician Dashboard.

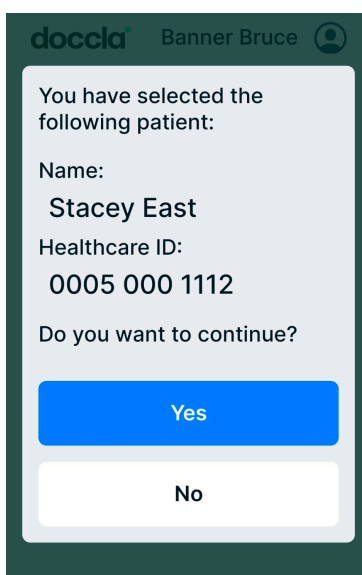

doccla Banner Bruce

You have selected the following patient:

Name:
Stacey East

Healthcare ID:
0005 000 1112

Do you want to continue?

Yes

No

Addendum

Reporting patient incident

Any serious incident that occurred in connection with the use of Doccla RPM must be reported to Doccla UK via support@doccla.com and the competent authority of the Member State where the user and/or patient is established.

Precautions

The information in this document is subject to change without notice.

Legal Notices

Doccla UK declares that the Doccla RPM software application is a medical device and is CE-marked in conformity with Regulation (EU) 2017/745 (MDR). The device complies with the applicable General Safety and Performance Requirements of the MDR. Conformity assessment was carried out by Kiwa Assurance B.V, NB no. 1912.

Manufacturer



Doccla UK Limited
184 Shepherds Bush Road
Hammersmith
London W6 7NL
United Kingdom



MedEnvoy Global BV
Prinses Margrietplantsoen 33, suite 123
The Hague
2595 AM
Netherlands



Date: 2026-06
End of life: 2028-06

Manufacturer responsibility

The manufacturer is only responsible for the software of Doccla RPM itself and therefore does not take responsibility for the functionality of the measurement devices used, the tablet or other devices that the Doccla RPM app can communicate with. The manufacturer ensures compliance with the listed Bluetooth devices in this document. No parts of this document may be reproduced or translated without the prior written permission of the manufacturer.

End of life (EOL)

The product's End of life (EOL) is 2 years after release. After EOL there is no product support. The date for EOL can be found on the product label.

The device (Doccla RPM software) works as an assisting tool and is not a replacement of treatments. It is recommended to avoid performing a measurement with moist or wet hands. A pen or other hardware tool can be used in case of touch screen challenge.